


Escape[®]

Nitinol Stone Retrieval Basket

The Escape Basket provides physicians with a versatile nitinol basket designed for use in the ureter and kidney. Should there be a need to reduce stone mass once the stone is captured, the Escape Basket with adaptor is designed to allow for simultaneous insertion and use of a 200 μ holmium laser fiber.

DESIGNED FOR USE WITH 200 μ HOLMIUM LASER FIBER

The Escape Basket can be used in combination with a 200 μ laser fiber in most semi-rigid and flexible ureteroscopes.

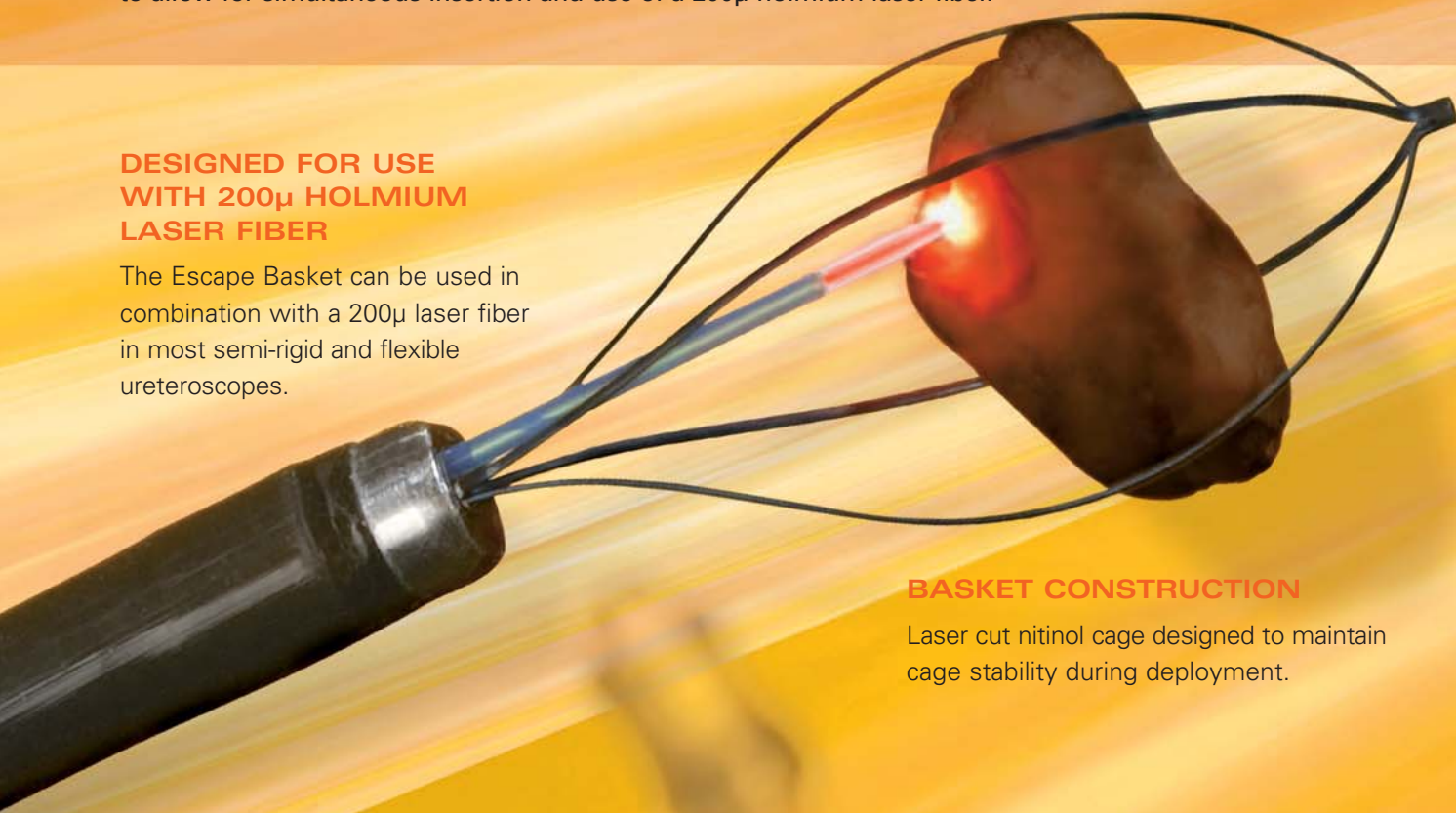

BASKET CONSTRUCTION

Laser cut nitinol cage designed to maintain cage stability during deployment.

**Boston
Scientific**

Delivering what's next.™

Escape®



Nitinol Stone Retrieval Basket

OUR BASKET DESIGN

- > In its fully open position, the 4 wire cage design extends 2 long legs for releasability.
- > In the standard 4-wire position the cage size is 11mm. In the fully extended position, the basket cage is designed to open to 15mm.

THE SIDE-BY-SIDE APPROACH

- > The Escape Basket is designed to allow a side-by-side approach, that is unique to the market, facilitating control and independent manipulation of both the basket and laser fiber.
- > With the basket and fiber side-by-side in the working channel of the scope, the laser fiber can be positioned between the basket wires for stone fragmentation and removal.
- > The Escape Basket is designed to stabilize the stone during lithotripsy, making it an effective anti-retropulsion device.*

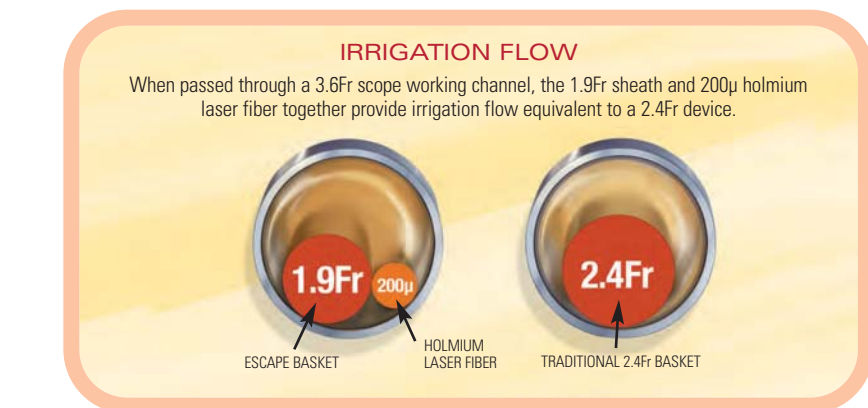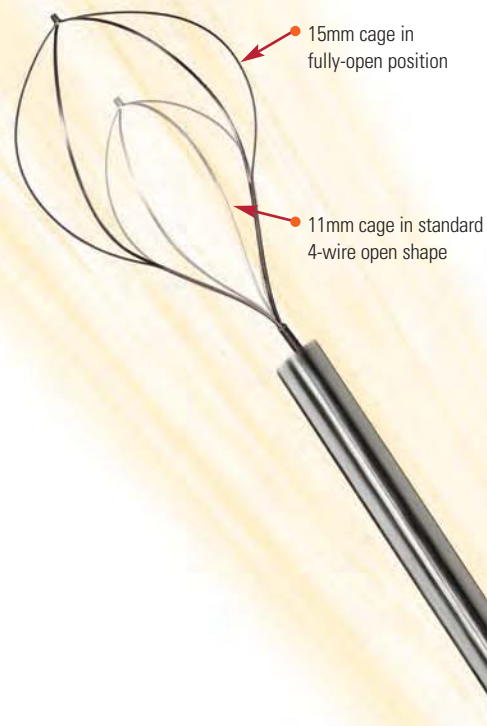

Ordering Information

ESCAPE NITINOL STONE RETRIEVAL BASKET			
Order Number	Description	Size	Qty.
M006390200	Escape Nitinol Stone Retrieval Basket	1.9Fr x 90cm	Each
M006390210	Escape Nitinol Stone Retrieval Basket	1.9Fr x 120cm	Each
OTHER RELATED PRODUCTS			
Order Number	Description		Qty.
M0068408921	200µ SlimLine EZ Single Use Holmium Laser Fiber		Box/3
M0068408400	200µ SlimLine Reusable Holmium Laser Fiber		Each

*Kesler et al., Use of the Escape® Nitinol Stone Retrieval Basket Facilitates Fragmentation and Extraction of Ureteral and Renal Calculi: A Pilot Study, *J of Endourology*, 2008;22(6):1213-7.

CAUTION: Federal Law (USA) restricts this device to sale by or on the order of a physician.
Refer to package insert provided with the product for complete Instructions for Use, Contraindications for Use, Potential Adverse Effects, Warnings and Precautions prior to using this product.

Boston Scientific

Delivering what's next.™

Boston Scientific Corporation
One Boston Scientific Place
Natick, MA 01760-1537
www.bostonscientific.com/urology

Ordering Information
1.888.272.1001

© 2009 Boston Scientific Corporation
or its affiliates. All rights reserved.

BVU2260 Rev.A 1M 12/09