


FilterWire EZ™

Embolic Protection System

Predictable Protection Made Easy*

Specifications

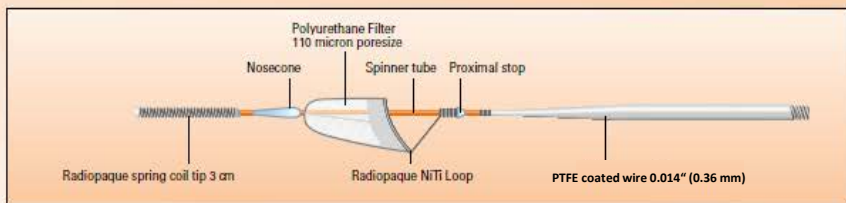

Ordering Information

Order Number	Description	Crossing profile	Retrieval profile	Vessel Size Coverage
H749 20100399 0	FilterWire EZ 190 cm	3.2F (1.07 mm)	4.3F (1.43 mm)	3.5 - 5.5 mm
H749 20100400 0	FilterWire EZ 300 cm	3.2F (1.07 mm)	4.3F (1.43 mm)	3.5 - 5.5 mm
H749 50100150 0	Bent Tip Retrieval Sheath EZ	4.0F (1.33 mm)	4.0F (1.33 mm)	N/A

* Data on file

**Boston
Scientific**
Delivering what's next.™

This material is not for use in the US or Japan.

All cited trademarks are the property of their respective owners. CAUTION: The law restricts these devices to sale by or on the order of a physician. Indications, contraindications, warnings and instructions for use can be found in the product labeling, included with each device. Information for the use of this device is available with applicable health authority product registrations.

Printed in Germany by medical/visio.

© 2007 Boston Scientific Corporation or its affiliates. All rights reserved.


IC-127801-AA

FilterWire EZ™

Embolic Protection System

**Boston
Scientific**

Predictable
Protection
Made Easy

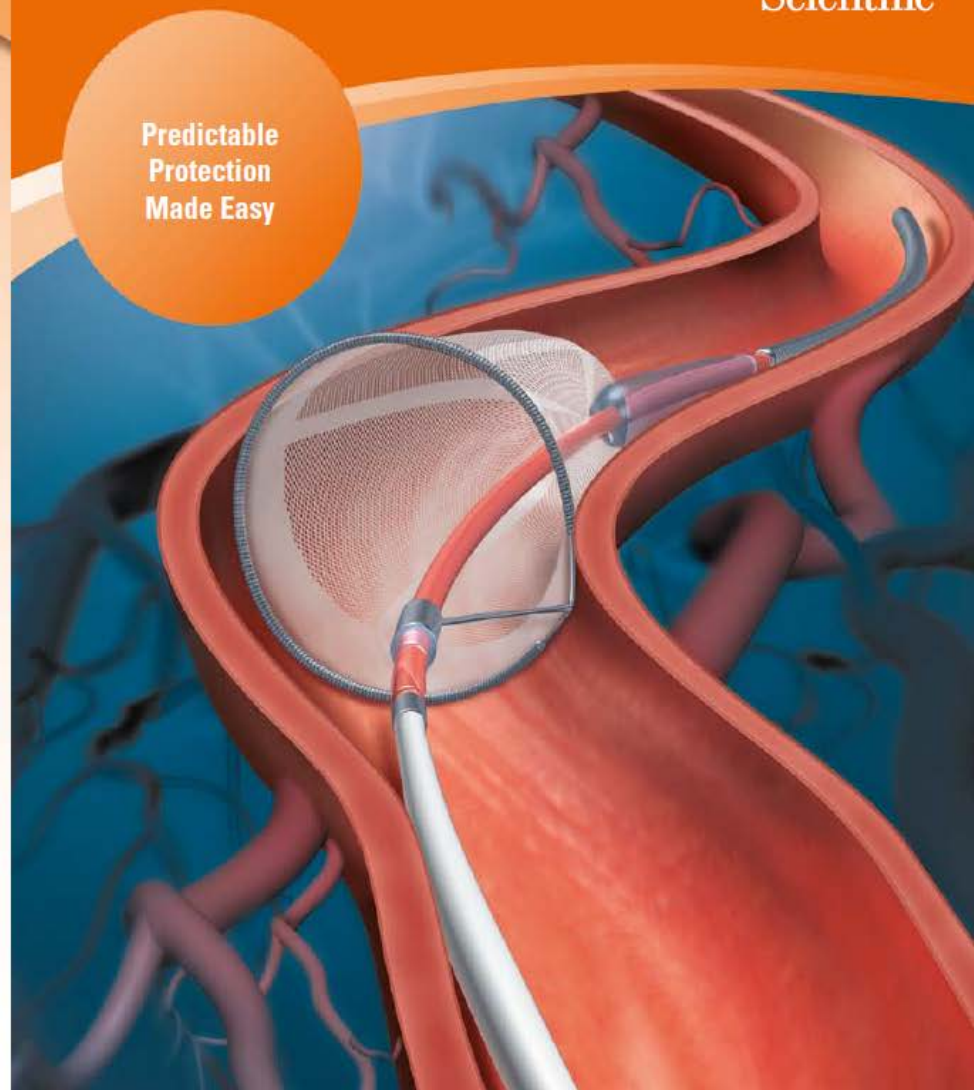


FilterWire EZ™

Embolic Protection System

Predictable Protection Made Easy*

FilterWire EZ™ adds predictability to your intervention by providing highly effective embolic protection without complicating your technique.

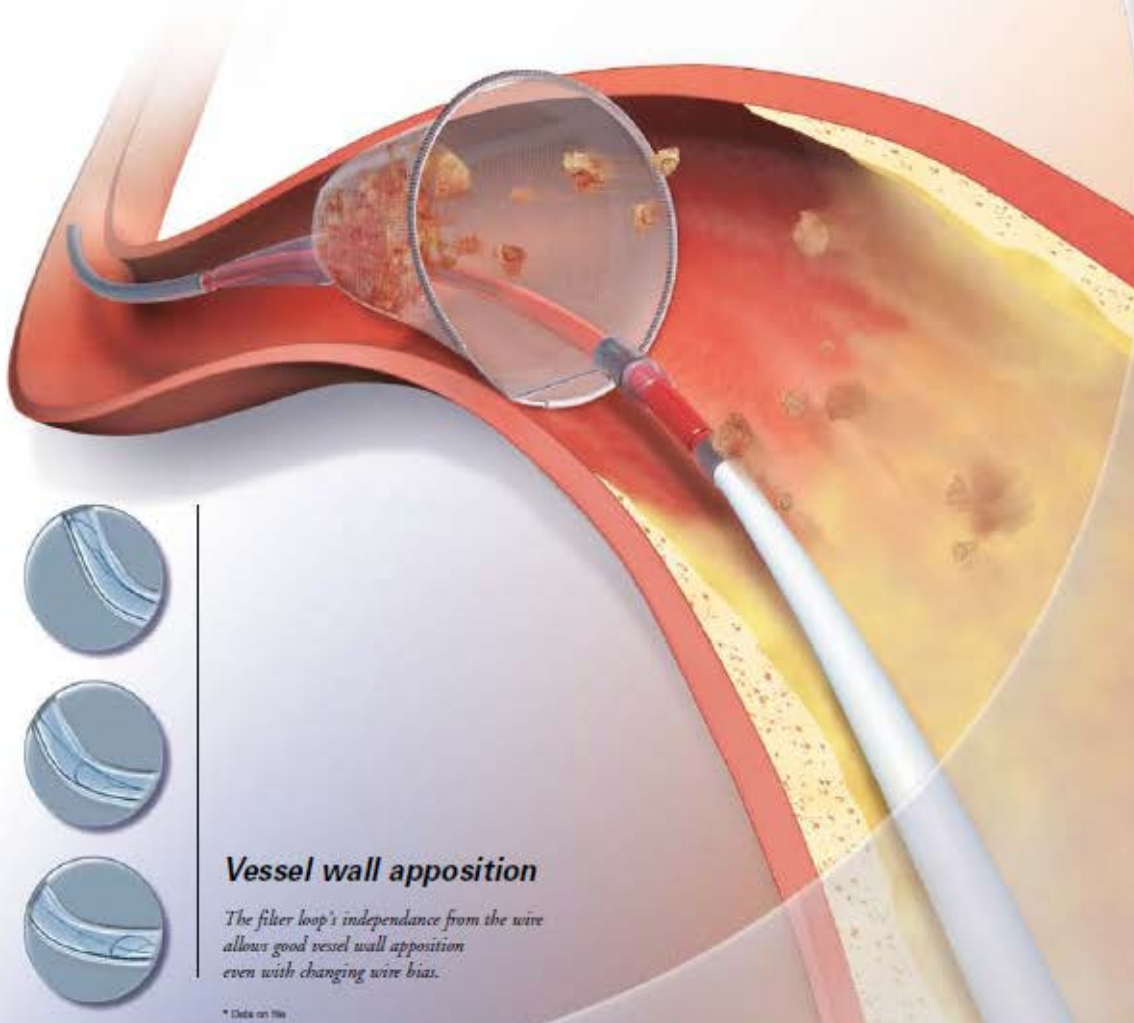

Vessel wall apposition

The filter loop's independence from the wire allows good vessel wall apposition even with changing wire bias.

* Click on file

Effective



The FilterWire EZ™ minimizes complications associated with percutaneous coronary interventions.



The suspended loop design allows the filter to be placed in curved vessels and still have good vessel wall apposition.



The conforming filter frame is designed to ensure reliable wall apposition in vessels of varying sizes and changes in the vessel.

Simple



A highly flexible tip with a small crossing profile and smooth transition makes the FilterWire EZ very deliverable.



The 3.5-5.5 mm artery size protection provides a broad range of vessel size coverage with one device.



The pre-loaded, peel-away delivery sheath is easy to prepare and provides rapid exchange convenience.